



MULTI SCIENCES

联科生物

Human IgA ELISA Kit

Catalog Number: EK174

Size: 48 Test, 96 Test, 2 × 96 Test, 5 × 96 Test, 10 × 96 Test

For the quantitative determination of human Immunoglobulin A (IgA) concentrations in cell culture supernates, serum and plasma.

This package insert must be read entirely before using this product. For proper performance, follow the protocol provided with each individual kit.

FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

MULTISCIENCES (LIANKE) BIOTECH, CO., LTD.

13F, Building 3, 3rd Phase, 108 Xiang Yuan Road, Gongshu Intellect
Information Industry Park, Hangzhou, Zhejiang Province, China.

www.multisciences.net

Tel: +86 057128828618-88662

Fax: +86-0571-28828618

E-mail: info.cn@liankebio.com

TABLE OF CONTENTS

ASSAY PROCEDURE SUMMARY 1

Introduction

Description 2
Principle of the Assay 2
Limitations of the Procedure 2

General Information

Materials Provided 3
Storage 3
Other Supplies Required 4
Precaution 4
Technical Hints 5









Assay Protocol

Sample Collection and Storage 5
Sample Preparation 6
Reagent Preparation 6, 7
Assay Procedure 7, 8

Analysis

Calculation of Results 8
Typical Data 8
Sensitivity 9
Precision 9
Recovery 9
Linearity 9
Calibration 10
Sample Values 10

ASSAY PROCEDURE SUMMARY

1. Prepare all reagents and standards as directed.

2. Add 50 μ l *Assay Buffer (1 \times)* and 50 μ l of 2-fold diluted *Standard* to Standard well in duplicate. Add 100 μ l *Assay Buffer (1 \times)* to Blank well in duplicate.

3. Add 90 μ l *Assay Buffer (1 \times)* and 10 μ l prediluted sample to the sample well (The dilution refers to the Sample Preparation on Page 6).

4. Add 50 μ l diluted *Detect Antibody* to each well. Step 2, 3 and 4 should be completed within 15 minutes.

5. Incubate for 3 hours at RT.

6. Aspirate and wash 6 times.

7. Add 100 μ l *Substrate Solution* to each well. Incubate for 5 - 30 minutes at RT. Protect from light.

8. Add 100 μ l *Stop Solution* to each well.

9. Read at 450 nm within 30 minutes. Correction 570 or 630 nm.

DESCRIPTION

Immunoglobulin A (IgA) is an antibody that plays a critical role in mucosal immunity. Its heavy chains are of the type α . IgA has two subclasses (IgA1 and IgA2) and can exist in a dimeric form called secretory IgA (sIgA). In its secretory form, IgA is the main immunoglobulin found in mucous secretions, it is also found in small amounts in blood. The secretory component of sIgA protects the immunoglobulin from being degraded by proteolytic enzymes, thus sIgA can survive in the harsh gastrointestinal tract environment and provide protection against microbes that multiply in body secretions. sIgA can also inhibit inflammatory effects of other immunoglobulins. IgA is a poor activator of the complement system, and opsonises only weakly.

PRINCIPLE OF THE ASSAY

This assay employs the quantitative sandwich enzyme immunoassay technique. A monoclonal antibody specific for human IgA has been pre-coated onto a microplate. Standard, samples and HRP-linked detect antibody specific for IgA are pipetted into the wells and IgA present is bound by the immobilized antibody and detect antibody following incubation. After washing, substrate solution is added to the wells and color develops in proportion to the amount of IgA bound in the initial step. The color development is stopped and the intensity of the color is measured.

LIMITATIONS OF THE PROCEDURE

- FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.
- Do not use expired kit or reagents.
- Do not use reagents from other lots or manufacturers. Do not prepare component by yourself.
- If concentration of assayed factor in samples is higher than the highest standard, dilute the serum/plasma samples with *Assay Buffer*, dilute the cell culture supernate samples with *cell culture medium*. Reanalyze these and multiply results by the appropriate dilution factor.
- Any variation in testing personnel, sample preparation, standard dilution, pipetting technique, washing techniques, incubation time, temperature, kit age and equipment can cause variation in results.
- This assay is designed to eliminate interference by factors present in biological samples. Until all factors have been tested in the ELISA immunoassay, the possibility of interference cannot be excluded.

MATERIALS PROVIDED (96 Test)

Unopened kit should be stored at 2 - 8°C.



- **IgA Microplate** (1 plate): 96-well polystyrene microplate (12 strips of 8 wells) coated with a monoclonal antibody against human IgA.
- **IgA Standard** (2 vials): Recombinant human IgA in a buffered protein base with preservatives; lyophilized.
- **IgA Detect Antibody** (1 vial, 80 µl): HRP-conjugate anti-human IgA detect antibody; 100× liquid.
- **Assay Buffer (10×)** (2 bottles, 10 ml): PBS with 0.5 % Tween-20 and 5 % BSA.
- **Substrate** (1 bottle, 15 ml): TMB (tetramethyl-benzidine).
- **Stop Solution** (1 bottle, 15 ml): 0.18 M sulfuric acid.
- **Washing Buffer (20×)** (1 bottle, 50 ml): PBS with 1 % Tween-20.
- **Plate Covers** (5 strips).

STORAGE

Store kit reagents between 2 and 8°C. Immediately after use remaining reagents should be returned to cold storage (2 to 8°C). Expiry of the kit and reagents is stated on labels.

Expiration date of the kit components can only be guaranteed if the components are stored properly, and if, in case of repeated use of one component, this reagent is not contaminated by the first handling.

Unopened kit		Store at 2 - 8°C (See expiration date on the label).
Opened/ Reconstituted Reagents	1× Washing Buffer 1× Assay Buffer Stop Solution Substrate TMB Detect Antibody	Up to 1 month at 2 - 8°C.
	Standard	Up to 1 month at ≤ -20°C in a manual defrost freezer. Discard after use.
	Microplate Wells	Up to 1 month at 2 - 8°C. Return unused strips to the foil pouch containing the desiccant pack, reseal along entire edge to maintain plate integrity.

Provided this is within the expiration date of the kit.

OTHER SUPPLIES REQUIRED

- **Microplate reader** capable of measuring absorbance at 450 nm, with correction wavelength set at 570 nm or 630 nm.
- **Pipettes and pipette tips.**
- 50 µl to 300 µl adjustable **multichannel micropipette** with disposable tips.
- Multichannel micropipette **reservoir**.
- **Beakers, flasks, cylinders** necessary for preparation of reagents.
- **Deionized or distilled water.**
- **Polypropylene** test tubes for dilution.

PRECAUTION

- All chemicals should be considered as potentially hazardous.
- We therefore recommend that this product is handled only by those persons who have been trained in laboratory techniques and that it is used in accordance with the principles of good laboratory practice. Wear suitable protective clothing such as laboratory overalls, safety glasses and gloves.
- Care should be taken to avoid contact with skin or eyes. In the case of contact with skin or eyes wash immediately with water. See material safety data sheet(s) and/or safety statement(s) for specific advice.
- The Stop Solution provided with this kit is an acid solution. Wear eyes, hand, face, and clothing protection when using this material.
- Reagents are intended for research use only and are not for use in diagnostic or therapeutic procedures.
- Do not mix or substitute reagents with those from other lots or other sources.
- Do not use kit reagents beyond expiration date on label.
- Do not expose kit reagents to strong light during storage and incubation.
- Do not eat or smoke in areas where kit reagents or samples are handled.
- Avoid contact of skin or mucous membranes with kit reagents or specimens.
- Rubber or disposable latex gloves should be worn while handling kit reagents or specimens.
- Avoid contact of substrate solution with oxidizing agents and metal.
- Avoid splashing or generation of aerosols.
- In order to avoid microbial contamination or cross- contamination of reagents or specimens which may invalidate the test use disposable pipette tips and/or pipettes.
- Use clean, dedicated reagent trays for dispensing the conjugate and substrate reagent.
- Exposure to acid inactivates the HRP and antibody conjugate.
- Glass-distilled water or deionized water must be used for reagent preparation.
- Substrate solution must be warmed to room temperature prior to use.
- Decontaminate and dispose specimens and all potentially contaminated materials as they could contain infectious agents. The preferred method of decontamination is autoclaving for a minimum of 1 hour at 121.5°C.
- Liquid wastes not containing acid and neutralized waste may be mixed with sodium hypochlorite in volumes such that the final mixture contains 1.0 % sodium hypochlorite. Allow 30 minutes for effective decontamination. Liquid waste containing acid must be neutralized prior to the addition of sodium hypochlorite.

TECHNICAL HINTS

- When mixing or reconstituting protein solutions, always avoid foaming.
- To avoid cross-contamination, change pipette tips between additions of each standard level, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.
- When using an automated plate washer, adding a 30 seconds soak period before washing step and/or rotating the plate between wash steps may improve assay precision.
- To ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary.
- Substrate Solution should remain colorless until added to the plate. Keep Substrate Solution protected from light. Substrate Solution should change from colorless to gradations of blue.
- Stop Solution should be added to the plate in the same order as the Substrate Solution.
- The color developed in the wells will turn from blue to yellow upon addition of the Stop Solution. Wells that are green in color indicate that the Stop Solution has not mixed thoroughly with the Substrate Solution.
- It is recommended that all samples and standards be assayed in duplicate.
- Take care not to scratch the inner surface of the microwells.

SAMPLE COLLECTION AND STORAGE

Cell Culture Supernates – Remove particulates by centrifugation at $300 \times g$ for 10 minutes and assay immediately or aliquot and store samples at $\leq -20^{\circ}\text{C}$.

Serum – Use a serum separator tube (SST) and allow samples to clot for 30 minutes before centrifugation for 10 minutes at $1,000 \times g$. Remove serum and assay immediately or aliquot and store samples at $\leq -20^{\circ}\text{C}$.

Plasma – Collect plasma using EDTA, citrate or heparin as anticoagulant. Centrifuge at $1,000 \times g$ within 30 minutes of collection. Assay immediately or aliquot and store samples at $\leq -20^{\circ}\text{C}$.

Other biological samples might be suitable for use in the assay. Cell culture supernates, serum and plasma were tested with this assay.

Note: Samples containing a visible precipitate must be clarified prior to use in the assay. Do not use grossly hemolyzed or lipemic specimens.

Samples should be aliquoted and must be stored frozen at -20°C to avoid loss of bioactive human IgA. If samples are to be run within 24 hours, they may be stored at 2 to 8°C .

Avoid repeated freeze-thaw cycles. Prior to assay, the frozen sample should be brought to room temperature slowly and mixed gently.

SAMPLE PREPARATION

Normal serum and plasma samples require a 10,000-fold dilution. A suggested 10,000-fold dilution is two step dilution: first, 10 μ l sample + 990 μ l Assay Buffer (1 \times); next, 10 μ l Mix + 990 μ l Assay Buffer (1 \times).

REAGENT PREPARATION

Bring all reagents and samples to room temperature before use.

If crystals form in the Buffer Concentrates, warm and gently stir them until completely dissolved.

Washing Buffer (1 \times)

Pour entire contents (50 ml) of the **Washing Buffer (20 \times)** into a clean 1000 ml graduated cylinder. Bring to final volume of 1000 ml with pure or deionized water.

Mix gently to avoid foaming.

Transfer to a clean wash bottle and store at 2 to 25 $^{\circ}$ C. Washing Buffer (1 \times) is stable for 30 days.

Assay Buffer (1 \times)

Pour the entire contents (20 ml) of the **Assay Buffer (10 \times)** into a clean 250 ml graduated cylinder. Bring to final volume of 200 ml with distilled water. Mix gently to avoid foaming.

Store at 2 to 8 $^{\circ}$ C. Assay Buffer (1 \times) is stable for 30 days.

Detect Antibody

Mix well prior to making dilutions.

Make a **1: 100** dilution of the concentrated **Detect Antibody** solution with Assay Buffer (1 \times) in a clean plastic tube as needed.

The diluted Detect Antibody should be used within 30 minutes after dilution.

Sample Dilution

If your samples have high IgA content, dilute serum/plasma samples with Assay Buffer (1 \times). For cell culture supernates, dilute with cell culture medium.

Human IgA Standard


Reconstitute **Human IgA Standard** by addition of distilled water. Reconstitution volume is stated on the label of the standard vial. Swirl or mix gently to insure complete and homogeneous solubilization (concentration of reconstituted standard = 800 ng/ml).

Allow the standard to reconstitute for 10 - 30 minutes. Mix well prior to making dilutions.

Use polypropylene tubes.

For serum/plasma samples, mixing *concentrated human IgA standard* (150 µl) with 150 µl of *Assay Buffer (1×)* creates the high standard (400 ng/ml). Pipette 150 µl of *Assay Buffer (1×)* into each tube. Use the high standard to produce a 1:1 dilution series (scheme below). Mix each tube thoroughly before the next transfer. *Assay Buffer (1×)* serves as the zero standard (0 pg/ml).

For cell culture supernates, mixing *concentrated human IgA standard* (150 µl) with 150 µl of cell culture medium creates the high standard (400 ng/ml). Pipette 150 µl of cell culture medium into each tube. Use the high standard to produce a 1:1 dilution series. Mix each tube thoroughly before the next transfer. Cell culture medium serves as the zero standard (0 pg/ml).

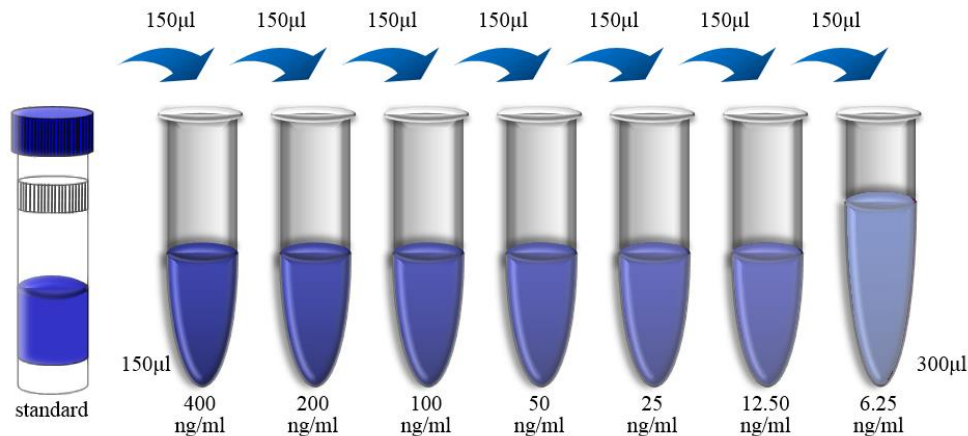

ASSAY PROCEDURE

Bring all reagents and samples to room temperature before use.

1. Prepare all reagents including microplate, samples, standards and working solution as described in the previous sections.
2. Remove excess microplate strips and return them to the foil pouch containing the desiccant pack, and reseal for further use.
3. Add 300 µl *Washing Buffer (1×)* per well, and allow it for about 30 seconds before aspiration. Soaking is highly recommended to obtain a good test performance. Empty wells and tap microwell strips on absorbent pad or paper towel to remove excess *Washing Buffer (1×)*. Use the microwell strips immediately after washing. **Do not allow wells to dry.**
4. Add 50 µl *Assay Buffer (1×)* and 50 µl of 2-fold diluted *Standard* to Standard well in duplicate. Add 100 µl of *Assay Buffer (1×)* to Blank well in duplicate.
5. Add 90 µl *Assay Buffer (1×)* and 10 µl of prediluted sample to the sample well (The dilution refers to the Sample Preparation on Page 6).
6. Add 50 µl of diluted *Detect Antibody* to each well. Ensure reagent addition in step 4, 5 and 6 is uninterrupted and completed within 15 minutes.
7. Cover with an adhesive strip. Incubate at room temperature (18 to 25°C) for 3 hours on a microplate shaker set at 300 rpm.
8. Aspirate each well and wash, repeating the process five times for a total six washes. Wash by filling each well with 300µl *Washing Buffer (1×)*. Complete removal of liquid at each step is essential to good performance. After the last wash, remove any remaining Wash Buffer by aspirating or decanting. Invert the plate and blot it against clean paper towels.

9. Add 100 µl of *Substrate Solution* to each well. Incubate for 5 - 30 minutes at room temperature. Protect from light.
10. Add 100 µl of *Stop Solution* to each well. The color will turn yellow. If the color in the well is green or if the color change does not appear uniform, gently tap the plate to ensure thorough mixing.
11. Measure the optical density value within 30 minutes by microplate reader set to 450 nm. If wavelength correction is available, set to 570 nm or 630 nm. If wavelength correction is not available, subtract readings at 570 nm or 630 nm from the readings at 450 nm. This subtraction will correct for optical imperfections in the plate. Reading directly at 450 nm without correction may generate higher concentration than true value.

CALCULATION OF RESULTS

Average the duplicate optical density readings for each standards and sample, then subtract the average optical density value of the zero standard.

Standard Concentration as horizontal axis, optical density (OD) Value as the vertical axis, regressing the data and create a standard curve using computer software. The data may be linearized by plotting the log of the IgA concentrations versus the log of the OD and the best fit line can be determined by regression analysis. This procedure will produce an adequate but less precise fit of the data.


Note: The finally concentration of top standard is 200 ng/ml. If instruction in this protocol have been followed samples have been diluted 1:9 (10 µl prediluted sample + 90 µl Assay Buffer), the concentration read from the standard curve must be multiplied by the dilution factor (×10).

If samples have been diluted following the instruction, the final dilution factor is 100,000. If sample have been diluted by other means, the concentration read from the standard curve must be multiplied by the appropriate dilution factor.

TYPICAL DATA

A standard curve must be run within each assay. This standard curve is provided for demonstration only.

ng/ml	O.D.	Average	Corrected
0.00	0.037	0.034	0.036
3.13	0.069	0.061	0.065
6.25	0.099	0.092	0.096
12.50	0.172	0.158	0.165
25.00	0.305	0.288	0.297
50.00	0.596	0.561	0.579
100.00	1.174	1.110	1.142
200.00	2.189	2.083	2.136

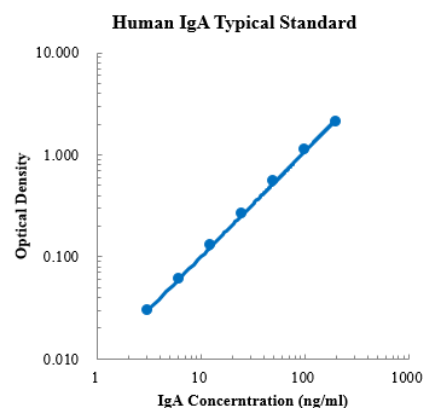

SENSITIVITY

The minimum detectable dose (MDD) of IgA is typically less than 0.35 ng/ml.

The MDD was determined by adding two standard deviations to the mean optical density value of ten zero standard replicates and calculating the corresponding concentration.

PRECISION

Intra-assay Precision (Precision within an assay)

Three serum-based and buffer-based samples of known concentration were tested twenty times on one plate to assess intra-assay precision.

Inter-assay Precision (Precision between assays)

Three serum-based and buffer-based samples of known concentration were tested in six separate assays to assess inter-assay precision.

	Intra-assay precision			Inter-assay precision		
Sample	1	2	3	1	2	3
n	20	20	20	6	6	6
Mean (ng/ml)	11.6	34.9	106.5	13.4	36.6	111.5
Standard deviation	0.6	1.2	5.0	0.3	0.3	7.6
CV (%)	4.8	3.4	4.7	2.4	0.7	6.8

RECOVERY

The spike recovery was evaluated by spiking 3 levels of human IgA into five health human serum samples. The un-spiked serum was used as blank in these experiments.

The recovery ranged from 82 % to 121 % with an overall mean recovery of 103 %.

LINEARITY

To assess the linearity of the assay, five samples were spiked with high concentration of IgA in human serum and diluted with Assay Buffer to produce samples with values within the dynamic range of the assay.

	Average (%)	Range (%)
1:2	96	90 - 108
1:4	99	88 - 109
1:8	97	91 - 106
1:16	99	95 - 108

CALIBRATION

This immunoassay is calibrated against a highly purified recombinant human IgA produced at MultiSciences.

SAMPLE VALUES

Serum/Plasma – Thirty samples from apparently healthy volunteers were evaluated for the presence of IgA in this assay. No medical histories were available for the donors used in this study.

Sample Matrix	Number of Samples Evaluated	Range (ug/ml)	Detectable (%)	Mean of Detectable (ug/ml)
Serum	30	289.3 - 2918.7	100	1623.1

Note: The sample range is non-physiological range. The sample range of healthy human will difference according to geographical, ethic, sample preparation, and testing personnel, equipment varies. The above information is only reference.

PLATE LAYOUT

